

Assembly Instructions

Airway Management Trainer P73 (1026007, 1026076)

1. Peel off the clear plastic film covering the suction base.



2. Gently pull back the neck skin of the model to expose the suction base lever.



3. Place the Airway Management Trainer on a clean, flat surface. Firmly press down the suction base and pull the lever downward to activate suction and lock the head securely in place.



4. Carefully pull back the neck skin to reveal the red circular connector inside the neck.



5. Insert the main bronchi tube into the connector and rotate it clockwise to lock it firmly.




6. Connect two of the green inflatable bags to the left and right main bronchi to simulate the lungs. Attach the third inflatable bag to the gastric port to simulate the stomach.



* Scan the barcode on the Smarter Simulation label to register your Airway Management Trainer and access the digital documents, spare parts and consumables information, product videos, and more.

* Please keep the label safe for any after-sales service.

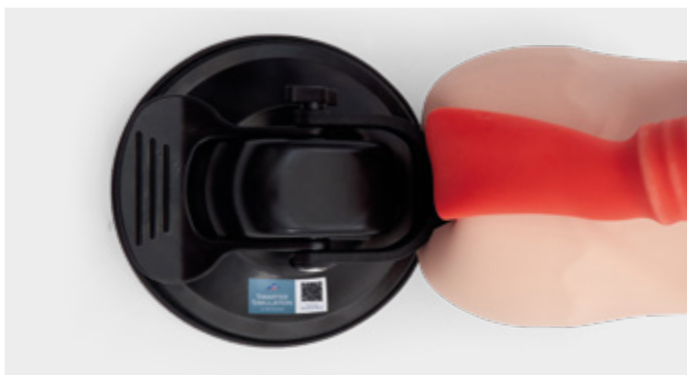

4th Order Segmental Bronchi (1026077)

The 4th Order Segmental Bronchi is an optional add-on, available for purchase separately, designed to enhance the functionality of the Airway Management Trainer.

How to use it: Detach the main bronchi from the red circular connector and attach the 4th Order Segmental Bronchi in its place.



IF YOU HAVE ANY QUESTIONS OR NEED SUPPORT, PLEASE SCAN HERE



3B Scientific Group

Headquarters · Ludwig-Erhard-Strasse 20 · 20459 Hamburg · Germany · 3bscientific.com