



Breast Examination, Biopsy and Ablation Training Platform



Instruction Manual

● Spare Parts

Full Model (BR-B01)



Single Breast Model
(BR-B02, BR-B03, BR-B06, BR-B07)




A. Platform base		1
B. Left single breast		1


C. Right single breast		1
D. Single breast support base		1

● Replace the Spare Parts (BR-B01)



- ① Lift up the left/right single breast using fingers from the platform base.

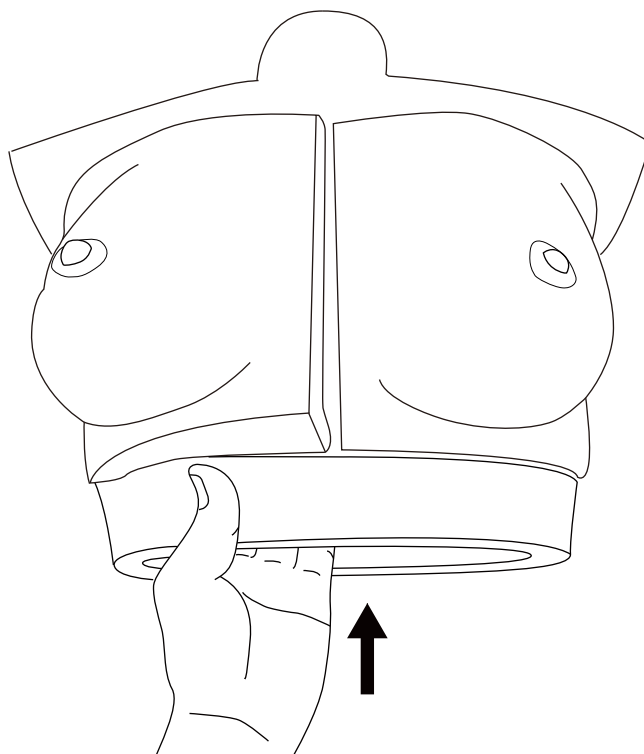

- ② Remove the left/right single breast.

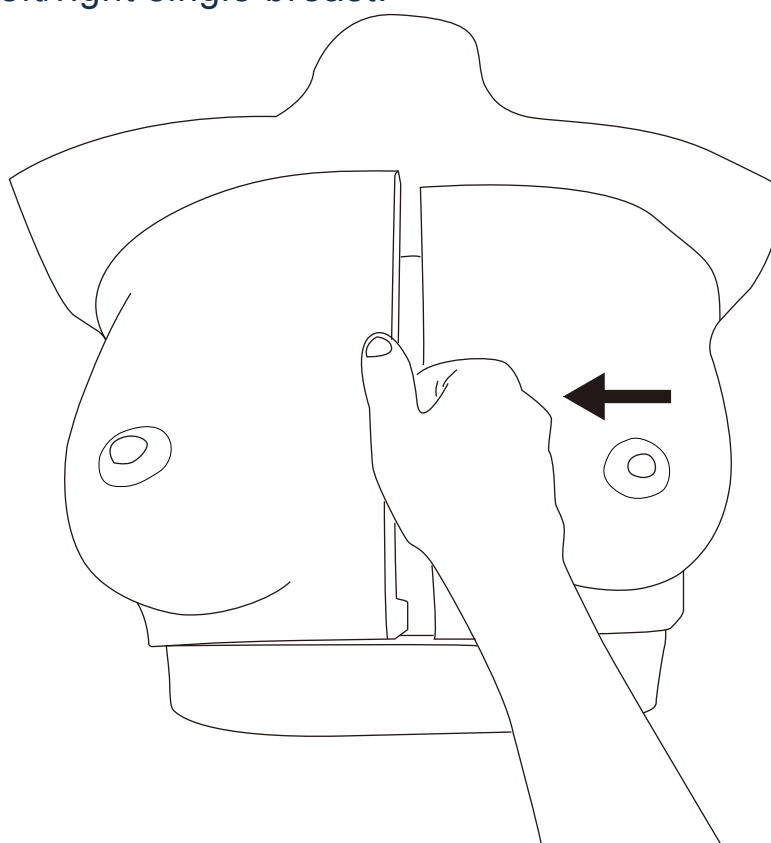

- ③ Replace and plug the left/right single breast on the frame.

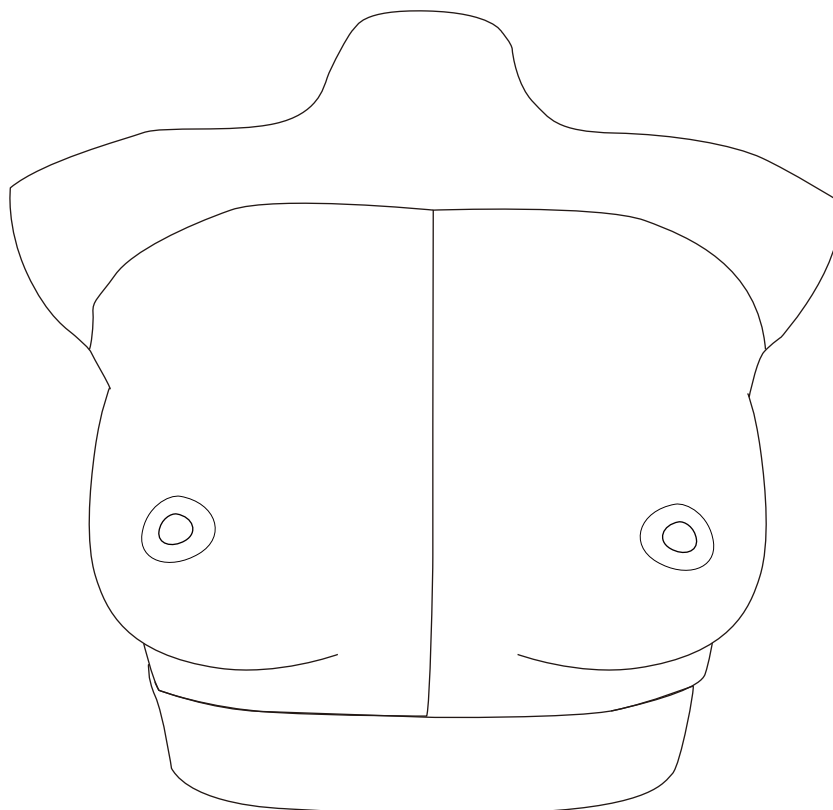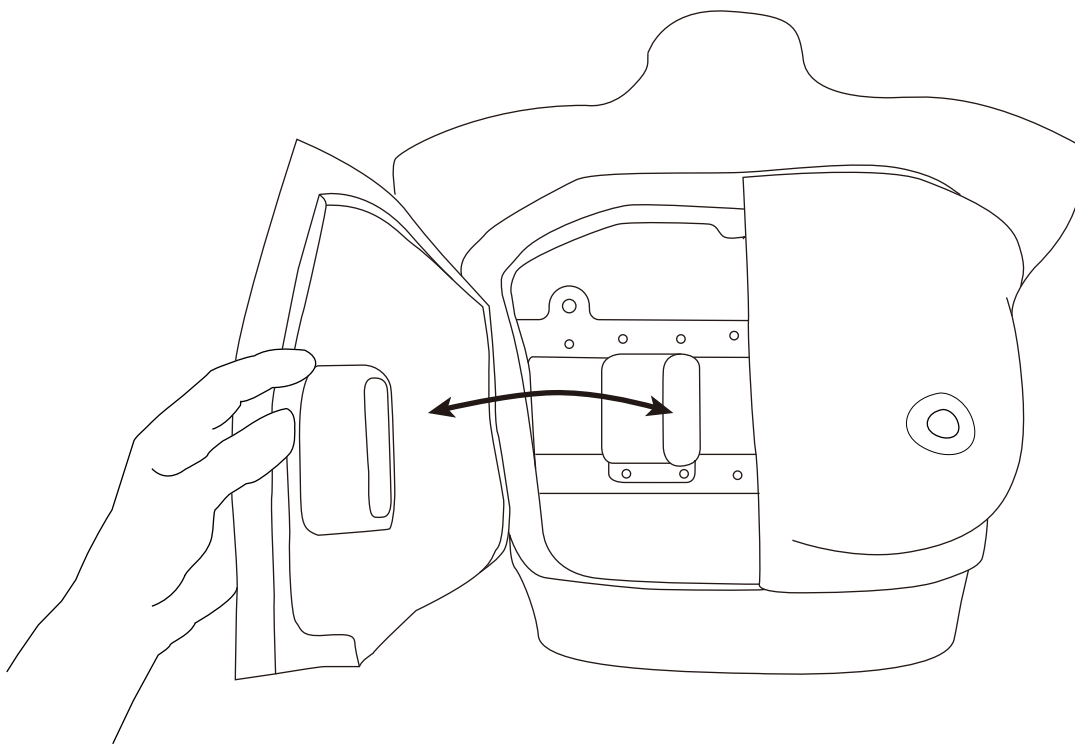

● Single Breast Model Setup

- ① Take out the single breast and the single breast support base (left or right, as applicable), and align them with the slots to assemble.



● Product Maintenance Guidelines

1. After each use (especially if ultrasound conductive gel has been applied), clean thoroughly and air dry at room temperature. Once dry, apply a thin layer of baby powder to maintain dryness.
2. Use a gentle, grease-removing detergent with simple ingredients, such as plain soap or mild dishwashing liquid, for cleaning.
3. Store the product in a cool, shaded indoor location.
4. The full model includes a protective portable case for transport. For long-distance shipping, consider adding an extra layer of bubble wrap around the platform base for additional protection.
5. The single breast model comes with cushioning materials inside the carton. Retain these materials and use them for protection during transport.



7F-7, No. 99, Sec. 1, Xintai 5th Rd.

Xizhi Dist., New Taipei City

221416, Taiwan

Tel: +886-2-2697-3933

E-mail: service@surgmate.com