


PRESTAN[®]

Ultralite[®]

with CPR Feedback

-User Guide-

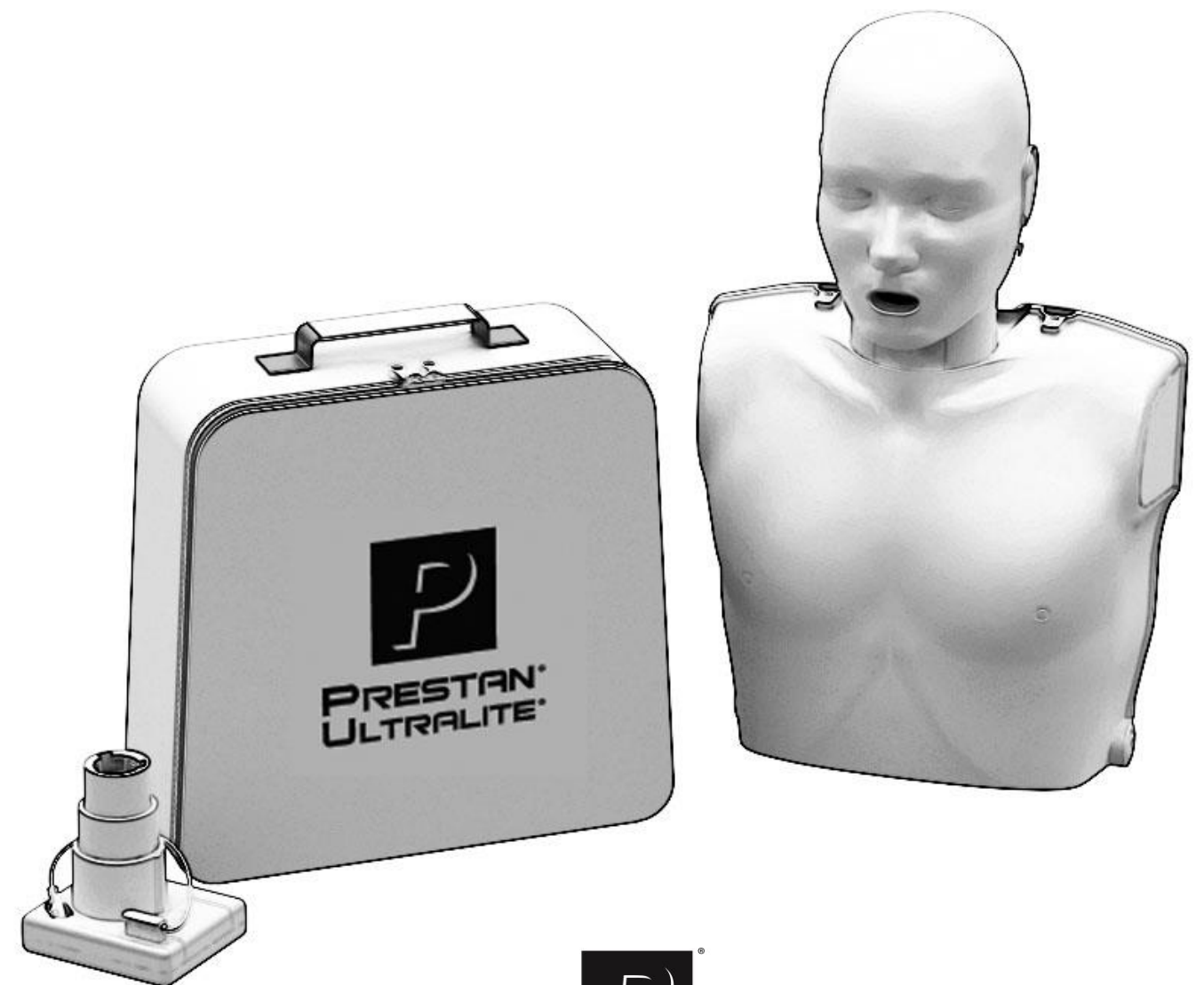



PRESTAN PRODUCTS THREE YEAR LIMITED WARRANTY

1. **WARRANTY COVERAGE:** Prestan offers a LIMITED WARRANTY against any defects in workmanship or materials under normal use with respect to its PRESTAN TRAINING MANIKIN.
2. **WARRANTY PERIOD:** Prestan warrants the PRESTAN TRAINING MANIKIN for three (3) years from the date of purchase. Any component of the PRESTAN TRAINING MANIKIN which proves defective within thirty-six (36) months from the date of purchase will be repaired or replaced, at Prestan's option, free of charge to the owner. The owner is responsible for shipping and handling fees to the place which Prestan designates.
3. **WARRANTOR'S RESPONSIBILITIES:** Prestan will repair or replace (at Prestan's option) a defective warranted component subject to the provisions of Section 2 without charge except for transportation costs. OWNER'S SOLE AND EXCLUSIVE REMEDY AGAINST PRESTAN SHALL BE FOR THE REPAIR OR REPLACEMENT OF DEFECTIVE GOODS AS PROVIDED HEREIN. Repair or replacement of any component of any PRESTAN TRAINING MANIKIN does not extend this limited warranty or begin a new limited warranty period.
4. **CIRCUMSTANCE NEGATING WARRANTY AND OTHER LIMITATIONS:**
 - A. Prestan reserves the right to refuse warranty service or replacement when product failure is due to misuse, mishandling, misapplication, accident, neglect, improper installation, modification, use of unauthorized parts or attachments, adjustments or repair performed by anyone other than Prestan or its authorized agent, or any causes other than defective workmanship or material.
 - B. This limited warranty does not cover abrasion, abnormal use, damage caused by mishandling, neglect or unauthorized modification. This warranty is void if the PRESTAN TRAINING MANIKIN is structurally altered.
 - C. Owner must provide proof of the date of purchase, place of purchase and the purchase price in order to obtain service under this limited warranty.
 - D. Prestan makes no warranty other than as contained herein. No agent, representative, dealer or employee of Prestan has the authority to alter the obligations or limitations of this limited warranty.
5. **LIMITATIONS OF IMPLIED WARRANTIES AND DAMAGES:**
 - A. Any and all implied warranties on the PRESTAN TRAINING MANIKIN, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, SHALL NOT EXCEED IN DURATION THE RESPECTIVE TERMS OF THIS LIMITED WARRANTY WHICH BEGINS WITH THE DATE OF PURCHASE. Some states do not allow limitations on how long an implied warranty lasts, so the above limitation may not apply to you.
 - B. PRESTAN SHALL NOT BE LIABLE FOR INCIDENTAL OR CONSEQUENTIAL DAMAGES ARISING OUT OF THE USE OF THE PRESTAN TRAINING MANIKIN OR THE BREACH OF THIS OR ANY OTHER EXPRESSED OR IMPLIED WARRANTY. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation or exclusion may not apply to you. Any and all damages awarded for any breach of this limited warranty shall be limited to the amount of the purchase price.
6. **WARRANTY SERVICE:** To obtain warranty service or information, contact Prestan Products, 701 Beta Drive, Suite 3, Mayfield Village, Ohio 44143 — or call Prestan toll free at 1-800-434-7989
7. **RIGHTS UNDER STATE LAW:** This warranty gives you specific legal rights and you may have other rights which vary from state to state.



PRESTAN® Ultralite® Manikins with CPR Feedback

THANK YOU for buying the Prestan Ultralite Manikins with CPR Feedback, # PP-ULM-100M-MS

CONTENTS:

- One (1) Manikin Torso
- One (1) Manikin Head Face
- One (1) Manikin Head Back
- One (1) CPR Feedback Piston (2 AA Batteries required, not included)
- One (1) Nylon Carrying Case
- One (1) 10-pack Ultralite Face-Shield/Lung-Bags

CARE AND CLEANING:

- When not in use, the components should be stored in the carrying case and in a cool, dry place.
- Ultralite Manikins may be cleaned in a variety of ways. A non-abrasive cleaner is recommended:
 - They can be wiped down with a Clorox® or Lysol® disinfecting wipe.
 - Cleaned with a mixture of soap and water.
 - Sprayed or wiped with common household cleaners such as Formula 409®, Windex® or Vim®.

The face may also be sanitized with alcohol, Clorox® or Lysol® wipes. Between student uses (and prior to applying the face shield), the instructor (or student) may clean the faces with a sanitizing wipe as an extra measure of protection.

Students may be provided with an individual face shield (sold separately) which may be placed over the Face-Shield/Lung-Bag. This allows multiple students to use one Face-Shield/Lung-Bag during a session.

Order through your Prestan Authorized Distributor:

Ultralite Manikins Face Shields, 50-count

#PP-UFS-50

Ultralite Manikins Face-Shield/Lung-Bags, 10-count

#PP-ULB-10

Ultralite Manikins Face-Shield/Lung-Bags, 50-count

#PP-ULB-50

TROUBLESHOOTING:

- If you cannot get the chest to rise when you blow into the face shield lung bag:
 - Check to make sure that the head is tilted back to create an open airway.
 - Check the lung bag itself. Ensure that it is properly positioned over the compression plate, not twisted and is lying flat.
 - Check the face shield. Ensure that it is attached to the 'buttons' on the side of the manikin head and that it is flat against the face of the manikin.
 - Check the bag to insure that has not been torn by testing inflation outside of the manikin.
- If the chest does not compress and/or it seems the piston is binding:
 - Ensure proper hand placement - check horizontal hand placement and adjust to center in between the nipples.

FURTHER INFORMATION:


Additional help guides and warranty information can be found at www.prestanproducts.com. If you have any questions about the Prestan Manikin line, if you would like to order more manikins, or if you would like to order any accessories for your manikin, please contact your Prestan Authorized Distributor.



PRESTAN®

Ultralite® Manikins with CPR Feedback








ASSEMBLE

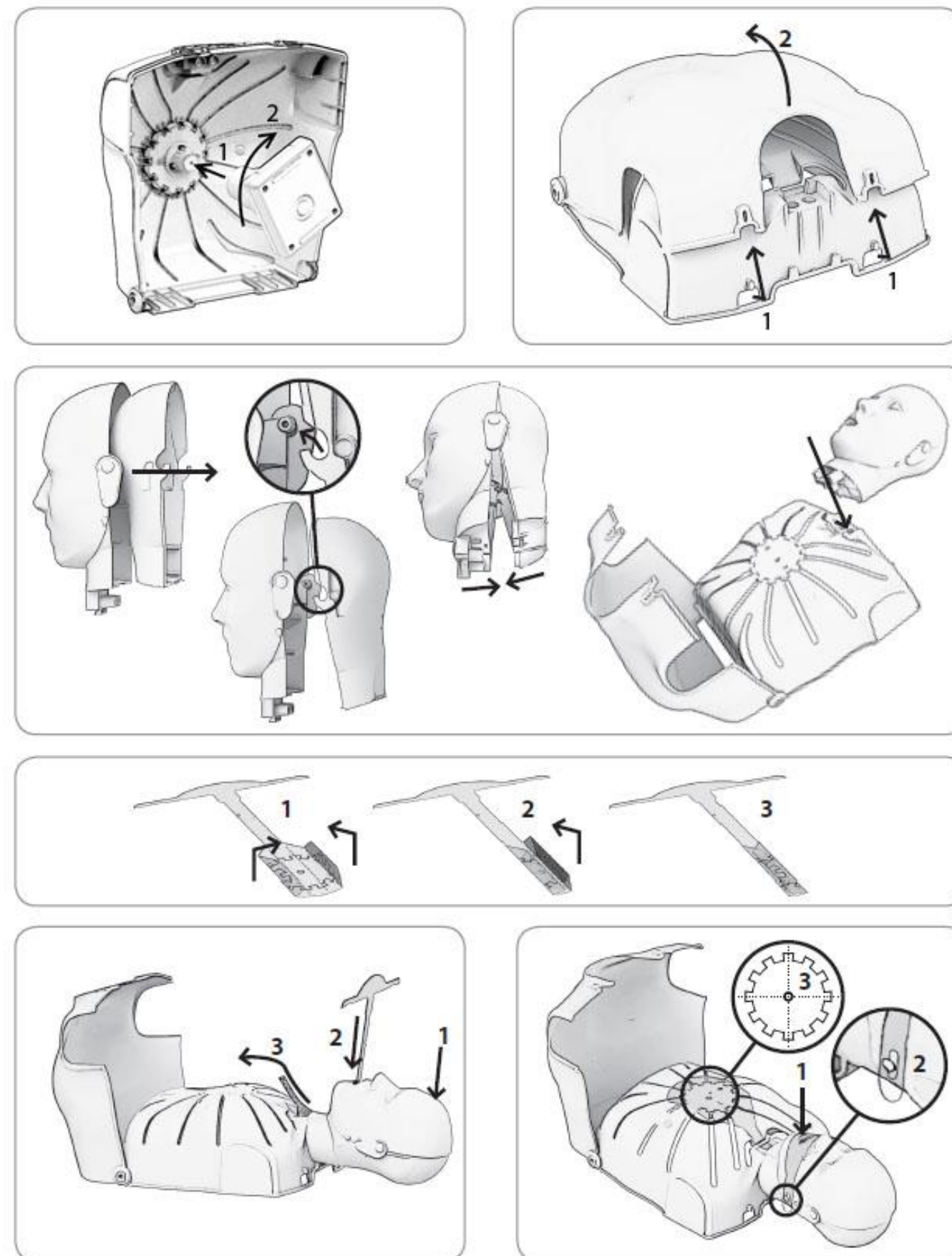

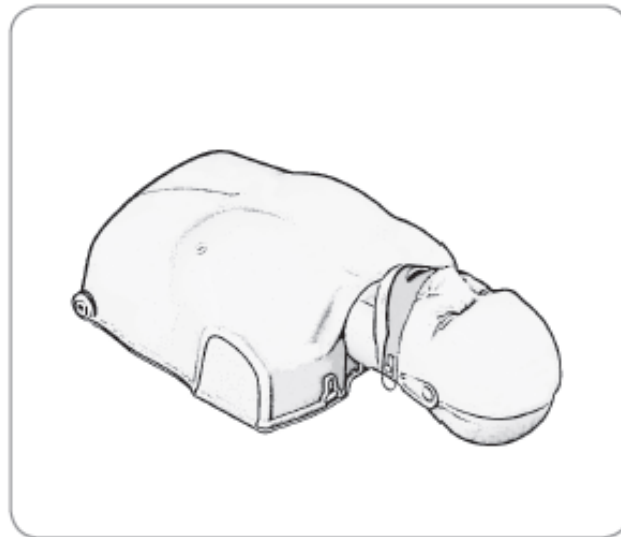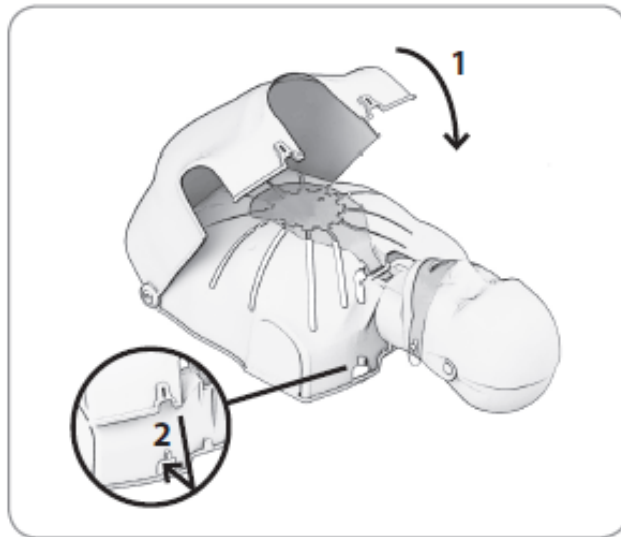

DISASSEMBLE

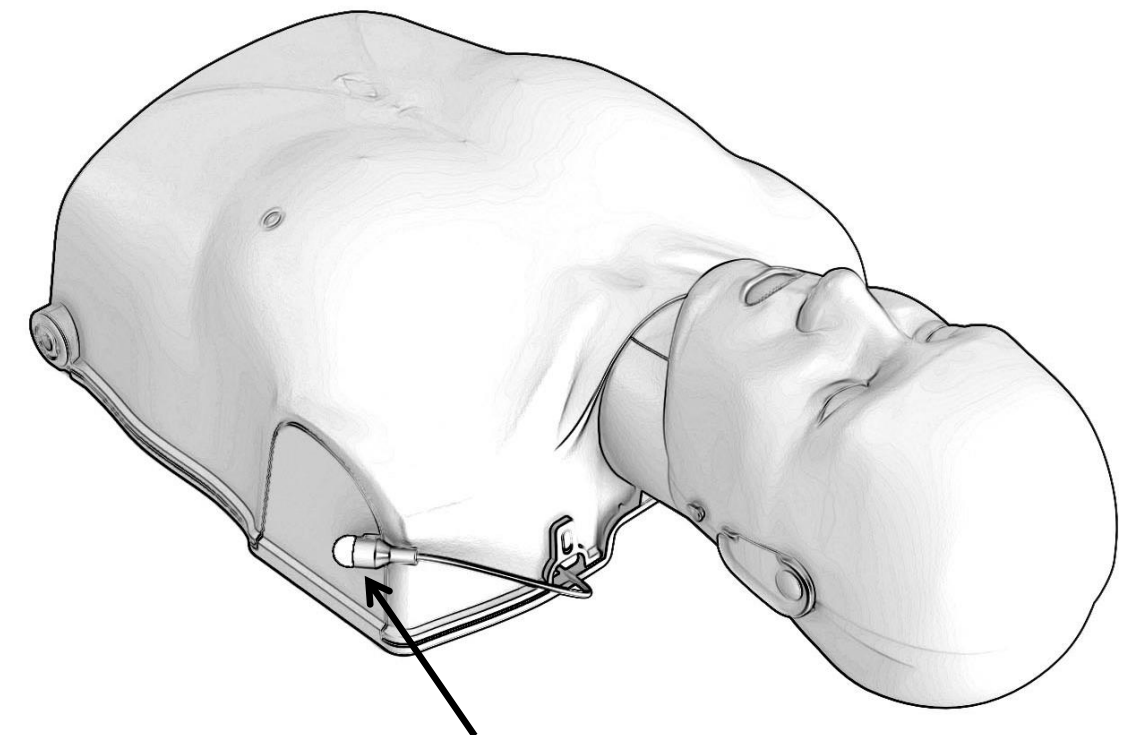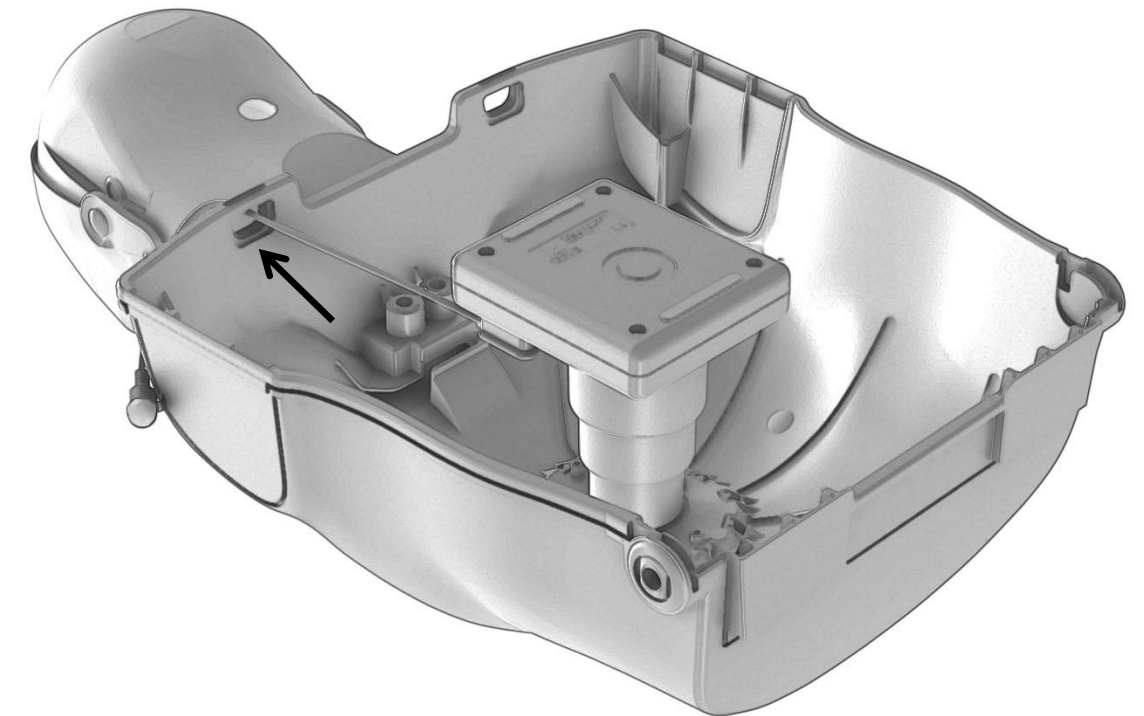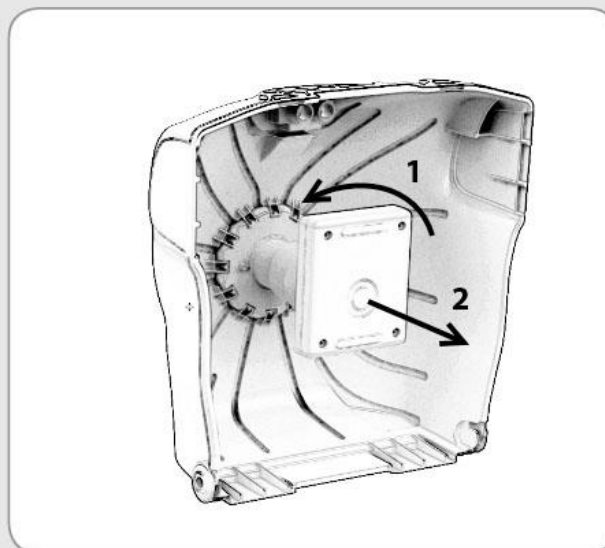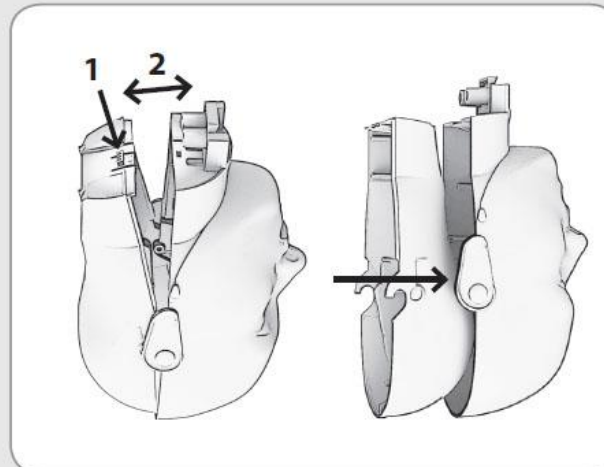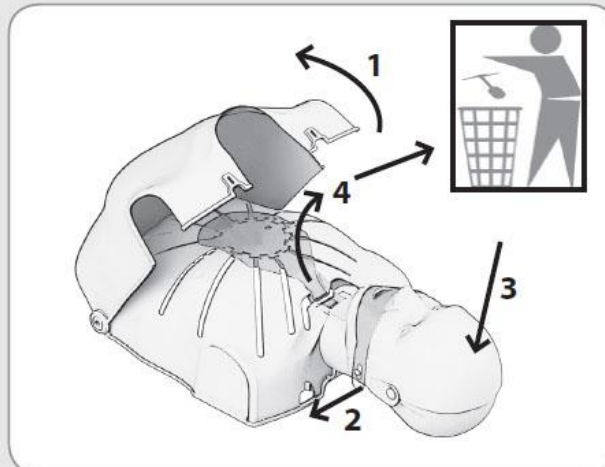

- Solid RED Light indicates under 100 CPM
- Solid GREEN Light indicates 100-120 CPM
(As recommended in 2015 AHA and ARC Guidelines)
- Flashing RED Light indicates more than 120 CPM

(CPM: Compressions per minute)



Prestan Professional Manikins

INFANT MANIKIN

USE AND CARE INSTRUCTION SHEET

THANK YOU for buying the **Prestan Professional Infant Manikin!**

Please read through these instructions and keep them with your manikin(s). For information on Prestan's Authorized Distributors, help guides, resources and more, visit the Prestan website at www.prestanproducts.com.



Prestan Professional Infant Manikin



Prestan Professional Infant Manikin 4-Pack

For future reference, make note below of the Prestan Authorized Distributor where you purchased your manikins. Call your Prestan Authorized Distributor for product questions, reordering parts/accessories, or warranty issues:

Prestan Authorized Distributor

Distributor Phone

Purchase Date



MODELS COVERED BY THIS INSTRUCTION SHEET

PRODUCT # (MODEL)	MANIKIN PACK	MANIKIN QTY./SIZE	INCLUDED CPR RATE MONITOR(S) *	INCLUDED CARRY CASE	INCLUDED FACE-SHIELD/ LUNG-BAGS
PP-IM-100M PP-IM-100M-MS PP-IM-100M-DS	SINGLE	1 INFANT	YES	SINGLE CASE	10-PACK
PP-IM-100 PP-IM-100-MS PP-IM-100-DS	SINGLE	1 INFANT	NO	SINGLE CASE	10-PACK
PP-IM-400M PP-IM-400M-MS PP-IM-400M-DS	4-PACK	4 INFANT	YES 1 PER MANIKIN	4-PACK CASE w/strap	50-PACK
PP-IM-400 PP-IM-400-MS PP-IM-400-DS	4-PACK	4 INFANT	NO	4-PACK CASE w/strap	50-PACK
PP-FM-300M PP-FM-300M-MS PP-FM-300M-DS	COLLECTION **	1 ADULT, 1 CHILD, 1 INFANT	YES 1 PER MANIKIN	COLLECTION CASE w/strap	10-PACK PER MANIKIN
PP-FM-500M PP-FM-500M-MS PP-FM-500M-DS	FAMILY PACK **	2 ADULTS, 1 CHILD, 2 INFANTS	YES 1 PER MANIKIN	FAMILY PACK CASE w/strap	10-PACK PER MANIKIN

REGARDING THE ABOVE CHART AND THE CONTENTS OF YOUR ORDER:

* The **CPR Rate Monitor** (*Rate Monitor, CPR Monitor, Monitor*) is a patented feedback feature for **Prestan Professional Manikins** that consists of an LED light bank that responds to compressions (see the section on the CPR Rate Monitor). The CPR Rate Monitor (*requires two (2) 'AA' batteries; not included*) is pre-installed in the diaper area of the Infant manikin. If you purchased the Prestan Professional Manikin without the CPR Rate Monitor and later decide you would like one, it may be ordered separately from your Prestan Authorized Distributor and easily installed. Installation instructions may be found in the Resources section of www.prestanproducts.com.

** **The Prestan Collection and Prestan Family Pack** also include **Prestan Professional Adult and Child Manikins**. These have separate Instruction Sheets that are included with the manikins or can be accessed at www.prestanproducts.com.

Available Skin Tones: Any of the above models may be available in one of three Skin Tones:

- A model number without a suffix signifies our original Light Skin Tone
- A model number with a “-MS” at the end (i.e. PP-AM-100M-MS) signifies a Medium Skin Tone
- A model number with a “-DS” at the end (i.e. PP-AM-400-DS) this signifies Dark Skin tone

If anything is missing, please contact the Prestan Authorized Distributor where you made your purchase.

CLEANING

- The face and the body skin of the Prestan® Professional Manikin is made out of the same material and can be cleaned in the same manner.
- If the material becomes dirty, clean with a mixture of soap and water. If the area is extremely dirty, common household cleaners such as Formula 409®, Windex® or Vim® may be used. Non-bleach wipes such as Clorox® or Lysol® disinfecting wipes are convenient and work well. It is recommended that the use of these cleaners be done cautiously with minimum amount of abrasion to insure the life of the skin.
- A melamine foam eraser such as Mr. Clean® Magic Eraser® or Scotch-Brite™ Easy Erasing Pad may be used cautiously for stains.
- Cleaners containing wax, oil, or citrus (d-Limonene) are not recommended.
- Any cleaner not listed/recommended above may be tested in a non-visible area of the manikin. Open the manikin's torso chest and test the cleaner on the underside. Monitor for 24 hours. Should there be an undesirable reaction, the cleaner should not be used on the manikin.
- The face skin may also be sanitized with alcohol wipes or soaked in a solution of 1/4 cup bleach mixed with one gallon of water for 10 minutes. Between student uses, the instructor (or student) may clean the faces off with an alcohol wipe as an extra measure of safety.
- Don't forget to wipe the upper chest foam and skin in the same way in order to sanitize the inside of the manikin as well.
- Replacement skins for the face and torso may be purchased from your Prestan Authorized Distributor.

STORAGE

- Store the Prestan Professional Manikin(s) in the nylon carry bag(s).
- For added protection, store the manikin inside the carry bag, and place the carry bag inside the original shipping box.
- Do not stack the manikins unless they are inside the original shipping box or other protective container.
- The nylon carry bag for the Prestan Professional Manikin(s) is not meant to store or transport belongings (i.e. books) in addition to the manikin(s) - please use the carry bag only for the intended manikin(s) to prevent damage to the bag.
- Store the manikin(s) in a cool, dry area, at temperatures between 50°-70° F and 50% relative humidity.
- For storage longer than 1 month, it is recommended that the batteries be removed from the CPR Rate Monitor (if the manikin has this feature).

CPR FEATURES

FACE-SHIELD/LUNG-BAG FOR THE VISIBLE CHEST RISE

- The Face-Shield/Lung-Bag (PP-ILB-10, PP-ILB-50) combines the Face Shield and the Lung Bag into a single, t-shaped airway for the manikin. Air goes in through the Face-Shield opening to inflate the Lung-Bag, which provides the visible chest rise in response to administered CPR breaths.
- Face-Shield/Lung-Bags are sized specifically to fit each Prestan Professional Manikin (Adult/Child/Infant) and are not interchangeable among the different manikin sizes.
- The Face-Shield portion covers the mouth of the manikin. It stays in place by hooking the ends onto buttons located at each side of the head under the ears.
- The Lung-Bag portion is placed on the manikin chest (under the torso skin) and fills with air from administered CPR breaths. This provides a visual chest rise for the student.
- Face-Shield/Lung-Bags are a consumable item and should be changed depending on how they are used:
 - If breaths are being administered by mouth, either directly or via rescue mask adaptor, the Face-Shield/Lung-Bag should be changed between classes due to accumulation of bacteria-causing moisture inside the neck and bag. This change is advised for both sanitary and performance purposes.
 - If breaths are being administered by a manual resuscitator (BVM/Ambu bag), there is no need to change the Face-Shield/Lung-Bag.
- When a manikin is being shared by more than one student, it is recommended that each student be provided with their own sanitary and disposable mouth protector:
 - If direct mouth contact is being used to give CPR breaths, students should be supplied with their own Face Shield (PP-IFS-50) (sold separately and disposable). The fitted Face Shield can be placed over the already-installed Face-Shield/Lung-Bag on the mouth of the manikin when it is that student's turn to give breaths. Generic, non-fitted face shields may also be used, but may slide around on the mouth and possibly cover the mouth opening if not held in place.
 - If rescue masks are being used to give CPR breaths, students should be supplied with their own Rescue Mask Adaptor (10076-PPA, 10076-PPA-50) (sold separately and disposable). The adaptors can be inserted into any standard mask when it is that student's turn to give breaths.

CPR CLICKER

- The Clicker is the round disk in the center of the manikin's back, under the infant skin and foam piece.
- The Clicker works with compressions to provide an audible "click" when compression depth reaches 2" (the compression is able to continue to a 2.4" depth).
- The Clicker sends a signal to the CPR Rate Monitor (if the manikin has this feature) to allow it to provide feedback on compression speed (see CPR Rate Monitor – Guide).

CPR RATE MONITOR - GUIDE

The CPR Rate Monitor (*requires two (2) 'AA' batteries; not included*) is a patented feedback feature that comes pre-installed in the diaper area of the Prestan Professional Child Manikin. If you have purchased a Prestan Professional Infant Manikin without the CPR Rate Monitor and later decide you want this feature – it may be purchased and installed separately by contacting your Prestan Authorized Distributor.

The Prestan Professional Training Manikin CPR Rate Monitor has a unique way of monitoring both the depth of compressions and the rate of compressions:

Press the button on the diaper to turn on the monitor. The smaller, green on-light will stay on until you press the button again to turn the monitor off. All four of the monitor lights will come on, indicating that the lights are working. (Figure 1)

If the lights do not light, make sure that the depth of compression is deep enough (you will hear a clicking sound). The lights will not work if the student is not compressing the manikin correctly.

Once a student begins chest compressions, from zero to less than 60 compressions per minute, the red light will be on. (Figure 2)

As the student increases the rate of compressions and reaches 60 compressions per minute, a yellow light will come on until the student reaches 80 compressions per minute. (Figure 3)

As the student increases the rate of compressions and reaches 80 compressions per minute, the first green light will come on until the student reaches 100 compressions per minute. (Figure 4)

When the student reaches 100 compressions per minute, a second green light comes on, indicating that the recommended rate has been reached. (Figure 5)

If the student slows down, one of the green lights will go out indicating that the student needs to increase the rate of compressions.

Should compressions exceed 120 per minute, a yellow light will flash rapidly indicating the recommended cap has been reached. (Figure 6)

Always try to maintain the recommended speed between 100 and 120 compressions per minute as indicated by the 2 green lights.

If the student's depth of compression is not correct, at any time, all the lights will go off and the red light will blink.



Figure 1



Figure 2; <60 cpm



Figure 3; 60-80 cpm



Figure 4; 80-100 cpm



Figure 5; 100-120 cpm



Figure 6; >120 cpm

FACE-SHIELD/LUNG BAG INSTALLATION

1. Unlatch the front torso skin of the Prestan Professional Infant Manikin by unhooking skin from the 'tabs' located at the top of the shoulders and peeling down the skin to expose the chest plate (Figure 1).
2. Fold the Lung-Bag portion of the Face-Shield/Lung-Bag (Figures 2a-2d):
 - The Face-Shield/Lung-Bag has four 'dotted lines' printed on one side of the bag (along with an oval)
 - Using the four dotted lines, fold the Face-Shield/Lung-Bag into fourths, following the diagram below



Figure 1

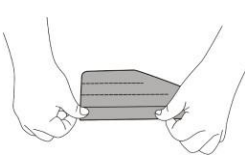


Figure 2a

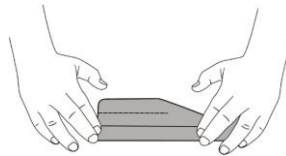


Figure 2b

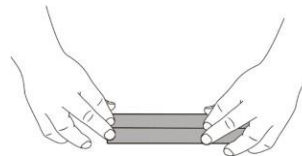


Figure 2c

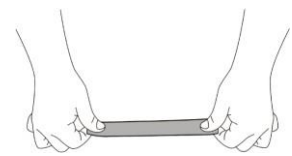


Figure 2d

3. Tilt back the head of the manikin to open the airway.
4. Insert the folded Lung-Bag into the mouth of the manikin (Figure 3).
5. Push the Lung-Bag into the mouth until the bag exits through the neck of the manikin (Figure 4).
6. Grasp the Lung-Bag and finish pulling it through the neck of the manikin (Figure 4).
7. Attach the Face-Shield tabs to the 'buttons' located on the side of the manikin head (Figure 5).
8. Lightly pull down on the Lung-Bag to help open the mouth hole on the Face-Shield.
9. Then position the Lung-Bag on top of the compression plate (the hard plastic circular piece found on top of the torso foam core) with the circular imprint on the Lung-Bag aligning with the compression plate (Figure 6).
10. With proper positioning, the neck of the Face-Shield/Lung-Bag will have slack. Make sure it is not twisted.
11. Roll the skin back up over the Lung-Bag being careful not to press down and "fold" the skin. The torso may then be closed by pushing on down on the shoulder tabs until they snap into place.

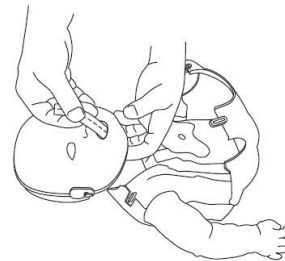


Figure 3

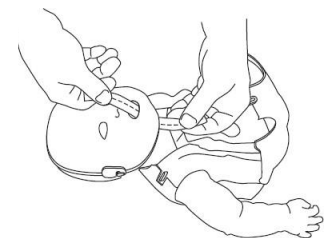


Figure 4

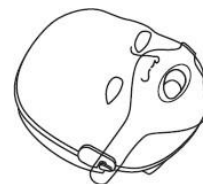


Figure 5

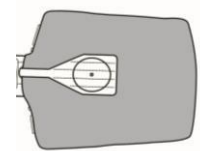


Figure 6

TROUBLESHOOTING GUIDE

SYMPTOM	POSSIBLE ISSUES OR COMMENTS	SOLUTIONS
CPR Rate Monitor		
Cannot find CPR Rate Monitor (Monitor) – no lights in the diaper area	You have purchased/received a Prestan Professional Infant Manikin without the CPR Rate Monitor	If you were expecting a monitored manikin, check with the distributor where you purchased the manikin to insure you received the item you paid for. See the chart on page 2 of this manual for the product numbers. Note: Infant CPR Rate Monitors may be purchased separately from your Prestan Authorized Distributor. They are shipped with a replacement diaper area cover for viewing the lights through it.
CPR Rate Monitor (Monitor) lights are not working properly <i>Please follow the troubleshooting tips at right in the order they are listed for best results</i>	- Chest compressions are not deep enough to signal Monitor - Batteries are not installed, need replacing or are not inserted correctly in the CPR Rate Monitor - Either the Clicker or the Monitor needs to be replaced	- Compress the chest until the click is audible from the Clicker at 2"-2.4" - Install or replace batteries into Monitor (back of manikin, left shoulder area) - requires 2-AA (1.5V) batteries - Exchange Monitors with working manikins (if available) to confirm whether the Monitor (RPP-IMON-1) needs to be replaced*
CPR Rate Monitor (Monitor) lights are all on	Batteries are running low	Replace batteries
CPR Rate Monitor (Monitor) red light is blinking	This is normal after replacing batteries or manikin is paused - may blink up to 16 times to indicate ready status	Commence or resume compressions - blinking will stop and feedback rate signals will begin
CPR Rate Monitor falls out when administering back blows	Monitor cover in diaper area is no longer fitting properly	Request new Infant Monitor Cover (10402)*

* Contact your Prestan Authorized Distributor to either purchase the replacement part or receive it at no charge if the manikin is within Prestan's 3-year warranty period and is a qualifying warranty repair.

TROUBLESHOOTING GUIDE

(Continued)

SYMPTOM	POSSIBLE ISSUES OR COMMENTS	SOLUTIONS
Torso and Face-Shield/Lung-Bag		
Cannot hear a click when doing compressions	- Hand position and/or compression angle may be incorrect - Compressions may not be deep enough to activate the clicker - Clicker may be broken	- Confirm correct hand placement and compress straight down from directly above the manikin - Compress between 2 and 2.4 inches depth - Replace Clicker in back of manikin (RPP-CCLICK-1)*
Chest will not rise when breaths are given	- Head is not tilted back enough to open airway - Face Shield is not attached - Lung Bag may be damaged or torn (Note: there are holes in the Infant Lung Bag necks by design to allow deflation) - Lung Bag is not lying flat or properly positioned over the compression plate and/or neck of bag is twisted - Lung Bag has accumulated moisture in the neck, which may cause the sides to stick together and prevent air passage	- Push on forehead and/or lift at chin to tilt head back enough to open airway - Attach Face Shield at “buttons” on the side of the manikin’s head below ears - Blow into the Face-Shield/Lung-Bag while it is outside the manikin to make sure it inflates properly – if not, discard and use another Face-Shield/Lung-Bag - Revisit instructions for installing a Face-Shield/Lung-Bag earlier in this manual - A new Face-Shield/Lung-Bag should be installed after each class, or more often if needed due to excessive moisture buildup
Manikin torso will not stay closed (when breaths are given)	Front torso skin may not be securely latched	Secure skin onto shoulder tabs

* Contact your Prestan Authorized Distributor to either purchase the replacement part or receive it at no charge if the manikin is within Prestan’s 3-year warranty period and is a qualifying warranty repair.

TROUBLESHOOTING GUIDE

(Continued)

SYMPTOM	POSSIBLE ISSUES OR COMMENTS	SOLUTIONS
Torso		
Rattling in chest	Compression Sleeve may be broken	Replace sleeve in chest with part RPP-ICOM-1 (compression assembly includes chest plate, spring, sleeve and clicker)*
Chest Plate is broken	Typically results from a combination of incorrect (not centered) hand placement, compression angle not straight down from directly above and/or overzealous compression force	Replace chest plate with RPP-ICOM-1 (compression assembly includes chest plate, spring, sleeve and clicker)*

* Contact your Prestan Authorized Distributor to either purchase the replacement part or receive it at no charge if the manikin is within Prestan's 3-year warranty period and is a qualifying warranty repair.

WARRANTY

Three Year Limited Warranty

1. **WARRANTY COVERAGE:** Prestan offers a LIMITED WARRANTY against any defects in workmanship or materials under normal use with respect to its PRESTAN TRAINING MANIKIN.
2. **WARRANTY PERIOD:** Prestan warrants the PRESTAN TRAINING MANIKIN for three (3) years from the date of purchase. Any component of the PRESTAN TRAINING MANIKIN which proves defective with thirty-six (36) months from the date of purchase will be repaired or replaced, at Prestan's option, free of charge to the owner. The owner is responsible for shipping and handling fees to the place which Prestan designates.
3. **WARRANTOR'S RESPONSIBILITIES:** Prestan will repair or replace (at Prestan's option) a defective warranted component subject to the provisions of Section 2 without charge except for transportation costs. OWNER'S SOLE AND EXCLUSIVE REMEDY AGAINST PRESTAN SHALL BE FOR THE REPAIR OR REPLACEMENT OF DEFECTIVE GOODS AS PROVIDED HEREIN. Repair or replacement of any component of any PRESTAN TRAINING MANIKIN does not extend this limited warranty or begin a new limited warranty period.
4. **CIRCUMSTANCE NEGATING WARRANTY AND OTHER LIMITATIONS:**
 - a. Prestan reserves the right to refuse warranty service or replacement when product failure is due to misuse, mishandling, misapplication, accident, neglect, improper installation, modification, use of unauthorized parts or attachments, adjustments or repair performed by anyone other than Prestan or its authorized agent, or any causes other than defective workmanship or material.
 - b. This limited warranty does not cover abrasion, abnormal use, damage caused by mishandling, neglect or unauthorized modification. This warranty is void if the PRESTAN TRAINING MANIKIN is structurally altered.
 - c. Owner must provide proof of the date of purchase, place of purchase and the purchase price in order to obtain service under this limited warranty.
 - d. Prestan makes no warranty other than as contained herein. No agent, representative, dealer or employee of Prestan has the authority to alter the obligations or limitations of this limited warranty.
5. **LIMITATIONS OF IMPLIED WARRANTIES AND DAMAGES:**
 - a. Any and all implied warranties on the PRESTAN TRAINING MANIKIN, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, SHALL NOT EXCEED IN DURATION THE RESPECTIVE TERMS OF THIS LIMITED WARRANTY WHICH BEGINS WITH THE DATE OF PURCHASE. Some states do not allow limitations on how long an implied warranty lasts, so the above limitation may not apply to you.
 - b. PRESTAN SHALL NOT BE LIABLE FOR INCIDENTAL OR CONSEQUENTIAL DAMAGES ARISING OUT OF THE USE OF THE PRESTAN TRAINING MANIKIN OR THE BREACH OF THIS OR ANY OTHER EXPRESSED OR IMPLIED WARRANTY. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation or exclusion may not apply to you. Any and all damages awarded for any breach of this limited warranty shall be limited to the amount of the purchase price.
6. **WARRANTY SERVICE.** To obtain warranty service contact your Prestan Authorized Distributor, or for information about Prestan, visit our website at www.prestanproducts.com.
7. **RIGHTS UNDER STATE LAW:** This warranty gives you specific legal rights and you may have other rights which vary from state to state.



Prestan Professional Manikins

ADULT MANIKIN

USE AND CARE INSTRUCTION SHEET

THANK YOU for buying the **Prestan Professional Adult Manikin!**

Please read through these instructions and keep them with your manikin(s). For information on Prestan's Authorized Distributors, help guides, resources and more, visit the Prestan website at www.prestanproducts.com.



Prestan Professional Adult Manikin

Prestan Professional Adult Manikin 4-Pack

For future reference, make note below of the Prestan Authorized Distributor where you purchased your manikins. Call your Prestan Authorized Distributor for product questions, reordering parts/accessories, or warranty issues:

Prestan Authorized Distributor

Distributor Phone

Purchase Date



Prestan Professional ADULT MANIKINS



MODELS COVERED BY THIS INSTRUCTION SHEET

PRODUCT # (MODEL)	MANIKIN PACK	MANIKIN QTY/SIZE	INCLUDED CPR RATE MONITOR(S) *	INCLUDED CARRY CASE	INCLUDED FACE-SHIELD/ LUNG-BAGS
PP-AM-100M PP-AM-100M-MS PP-AM-100M-DS	SINGLE	1 ADULT	YES	SINGLE CASE	10-PACK
PP-AM-100 PP-AM-100-MS PP-AM-100-DS	SINGLE	1 ADULT	NO	SINGLE CASE	10-PACK
PP-AM-400M PP-AM-400M-MS PP-AM-400M-DS	4-PACK	4 ADULTS	YES 1 PER MANIKIN	4-PACK CASE w/strap	50-PACK
PP-AM-400 PP-AM-400-MS PP-AM-400-DS	4-PACK	4 ADULTS	NO	4-PACK CASE w/strap	50-PACK
PP-FM-300M PP-FM-300M-MS PP-FM-300M-DS	COLLECTION **	1 ADULT, 1 CHILD, 1 INFANT	YES 1 PER MANIKIN	COLLECTION CASE w/strap	10-PACK PER MANIKIN
PP-FM-500M PP-FM-500M-MS PP-FM-500M-DS	FAMILY PACK **	2 ADULTS, 1 CHILD, 2 INFANTS	YES 1 PER MANIKIN	FAMILY PACK CASE w/strap	10-PACK PER MANIKIN

REGARDING THE ABOVE CHART AND THE CONTENTS OF YOUR ORDER:

* The **CPR Rate Monitor** (*Rate Monitor, CPR Monitor, Monitor*) is a patented feedback feature for **Prestan Professional Manikins** that consists of an LED light bank that responds to compressions (see the section on the CPR Rate Monitor). The CPR Rate Monitor (*requires two (2) 'AA' batteries; not included*) is pre-installed in the left shoulder of the adult manikin. If you purchased the Prestan Professional Manikin without the CPR Rate Monitor and later decide you would like one, it may be ordered separately from your Prestan Authorized Distributor and easily installed. Installation instructions may be found in the Resources section of www.prestanproducts.com.

** **The Prestan Collection** and **Prestan Family Pack** also include **Prestan Professional Infant and Child Manikins**. These have separate Instruction Sheets that are included with the manikins or can be accessed at www.prestanproducts.com.

Available Skin Tones: Any of the above models may be available in one of three Skin Tones:

- A model number without a suffix signifies our original Light Skin Tone
- A model number with a “-MS” at the end (i.e. PP-AM-100M-MS) signifies a Medium Skin Tone
- A model number with a “-DS” at the end (i.e. PP-AM-400-DS) this signifies Dark Skin tone

Black Felt Dots are included on a strip for optional use with Adult Manikins. These may be adhered to the underside of the manikin torso to help prevent sliding on slippery floor surfaces.

If anything is missing, please contact the Prestan Authorized Distributor where you made your purchase.

Prestan Professional ADULT MANIKINS



CLEANING

- The face and the body skin of the Prestan® Professional Manikin is made out of the same material and can be cleaned in the same manner.
- If the material becomes dirty, clean with a mixture of soap and water. If the area is extremely dirty, common household cleaners such as Formula 409®, Windex® or Vim® may be used. Non-bleach wipes such as Clorox® or Lysol® disinfecting wipes are convenient and work well. It is recommended that the use of these cleaners be done cautiously with minimum amount of abrasion to insure the life of the skin.
- A melamine foam eraser such as Mr. Clean® Magic Eraser® or Scotch-Brite™ Easy Erasing Pad may be used cautiously for stains.
- Cleaners containing wax, oil, or citrus (d-Limonene) are not recommended.
- Any cleaner not listed/recommended above may be tested in a non-visible area of the manikin. Open the manikin's torso chest and test the cleaner on the underside. Monitor for 24 hours. Should there be an undesirable reaction, the cleaner should not be used on the manikin.
- The face skin may also be sanitized with alcohol wipes or soaked in a solution of 1/4 cup bleach mixed with one gallon of water for 10 minutes. Between student uses, the instructor (or student) may clean the faces off with an alcohol wipe as an extra measure of safety.
- Don't forget to wipe the upper chest foam and skin in the same way in order to sanitize the inside of the manikin as well.
- Replacement skins for the face and torso may be purchased from your Prestan Authorized Distributor.

STORAGE

- Store the Prestan Professional Manikin(s) in the nylon carry bag(s).
- The manikin head should be removed for protection of the facial features, specifically the nose.
- For added protection, store the manikin inside the carry bag, and place the carry bag inside the original shipping box.
- Do not stack the manikins unless they are inside the original shipping box or other protective container.
- The nylon carry bag for the Prestan Professional Manikin(s) is not meant to store or transport belongings (i.e. books) in addition to the manikin(s) - please use the carry bag only for the intended manikin(s) to prevent damage to the bag.
- Store the manikin(s) in a cool, dry area, at temperatures between 50°-70° F and 50% relative humidity.
- For storage longer than 1 month, it is recommended that the batteries be removed from the CPR Rate Monitor (if the manikin has this feature).

CPR FEATURES

FACE-SHIELD/LUNG-BAG FOR THE VISIBLE CHEST RISE

- The Face-Shield/Lung-Bag (PP-ALB-10, PP-ALB-50) combines the Face Shield and the Lung Bag into a single, t-shaped airway for the manikin. Air goes in through the Face-Shield opening to inflate the Lung-Bag, which provides the visible chest rise in response to administered CPR breaths.
- Face-Shield/Lung-Bags are sized specifically to fit each Prestan Professional Manikin (Adult/Child/Infant) and are not interchangeable among the different manikin sizes.
- The Face-Shield portion covers the mouth of the manikin. It stays in place by hooking the ends onto buttons located at each side of the head under the ears.
- The Lung-Bag portion is placed on the manikin chest (under the torso skin) and fills with air from administered CPR breaths. This provides a visual chest rise for the student.
- Face-Shield/Lung-Bags are a consumable item and should be changed depending on how they are used:
 - If breaths are being administered by mouth, either directly or via rescue mask adaptor, the Face-Shield/Lung-Bag should be changed between classes due to accumulation of bacteria-causing moisture inside the neck and bag. This change is advised for both sanitary and performance purposes.
 - If breaths are being administered by a manual resuscitator (BVM/Ambu bag), there is no need to change the Face-Shield/Lung-Bag.
- When a manikin is being shared by more than one student, it is recommended that each student be provided with their own sanitary and disposable mouth protector:
 - If direct mouth contact is being used to give CPR breaths, students should be supplied with their own Face Shield (PP-AFS-50) (sold separately and disposable). The fitted Face Shield can be placed over the already-installed Face-Shield/Lung-Bag on the mouth of the manikin when it is that student's turn to give breaths. Generic, non-fitted face shields may also be used, but may slide around on the mouth and possibly cover the mouth opening if not held in place.
 - If rescue masks are being used to give CPR breaths, students should be supplied with their own Rescue Mask Adaptor (10076-PPA, 10076-PPA-50) (sold separately and disposable). The adaptors can be inserted into any standard mask when it is that student's turn to give breaths.

CPR CLICKER

- The Clicker is the round disk in the center of the manikin's back.
- The Clicker works with compressions to provide an audible "click" when compression depth reaches 2" (the compression is able to continue to a 2.4" depth).
- The Clicker sends a signal to the CPR Rate Monitor (if the manikin has this feature) to allow it to provide feedback on compression speed (see CPR Rate Monitor – Guide)

Prestan Professional ADULT MANIKINS



CPR RATE MONITOR - GUIDE

The CPR Rate Monitor (*requires two (2) 'AA' batteries; not included*) is a patented feedback feature that comes pre-installed in the shoulder area of the Prestan Professional Adult Manikin. If you have purchased a Prestan Professional Adult Manikin without the CPR Rate Monitor and later decide you want this feature – it may be purchased and installed separately by contacting your Prestan Authorized Distributor.

The Prestan Professional Training Manikin CPR Rate Monitor has a unique way of monitoring both the depth of compressions and the rate of compressions:

At the first compression all four of the lights in the upper left shoulder will come on, indicating that the lights are working. (Figure 1)

If the lights do not light, make sure that the depth of compression is deep enough (you will hear a clicking sound). The lights will not work if the student is not compressing the manikin correctly.

Once a student begins chest compressions, the yellow light and the green lights will go out and, from zero to less than 60 compressions per minute, the red light will be on. (Figure 2)

As the student increases the rate of compressions and reaches 60 compressions per minute, a yellow light will come on until the student reaches 80 compressions per minute. (Figure 3)

As the student increases the rate of compressions and reaches 80 compressions per minute, the first green light will come on until the student reaches 100 compressions per minute. (Figure 4)

When the student reaches 100 compressions per minute, a second green light comes on, indicating that the recommended rate has been reached. (Figure 5)

If the student slows down, one of the green lights will go out indicating that the student needs to increase the rate of compressions.

Should compressions exceed 120 per minute, a yellow light will flash rapidly indicating the recommended cap has been reached. (Figure 6)

Always try to maintain the recommended speed between 100 and 120 compressions per minute as indicated by the 2 green lights.

If the student's depth of compression is not correct, at any time, all the lights will go off and the red light will blink.



Figure 1

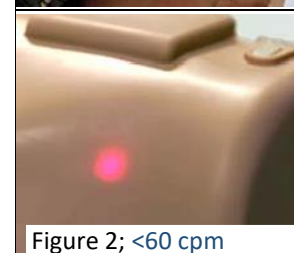


Figure 2; <60 cpm



Figure 3; 60-80 cpm



Figure 4; 80-100 cpm



Figure 5; 100-120 cpm

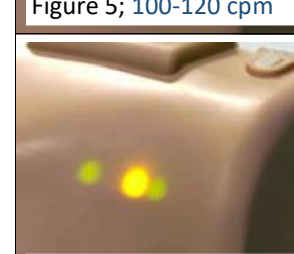


Figure 6; >120 cpm

FACE-SHIELD/LUNG-BAG INSTALLATION

1. Unlatch the torso of the Prestan Professional Manikin by popping the 'tabs' located at the top of the shoulders (Figure 1). Carefully lift the torso frame open by the neck.
Note: Do not pull on tabs (past popping the latch). Doing so may break the tabs or stretch out the frame over time, impairing the latching function.

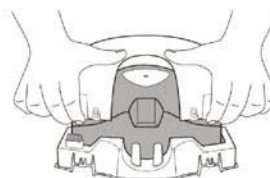


Figure 1



Figure 2

2. Attach the head to the torso by placing the neck of the head onto the two spindles on the neck portion of the torso (Figure 2).
3. Fold the Lung-Bag portion of the Face-Shield/Lung-Bag (Figures 3a-3d):
 - The Face-Shield/Lung-Bag has four 'dotted lines' printed on one side of the bag (along with an oval)
 - Using the four dotted lines, fold the Face-Shield/Lung-Bag into fourths, following the diagram below

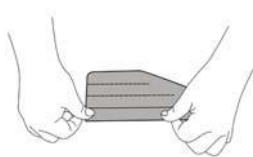


Figure 3a

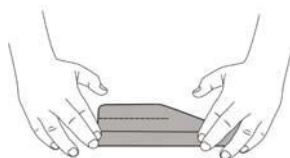


Figure 3b

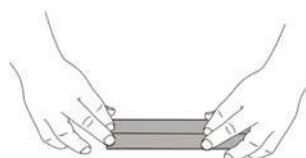


Figure 3c

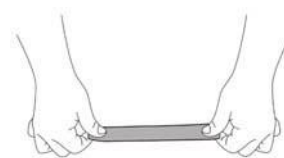


Figure 3d

4. Tilt back the head of the manikin (Figure 4).
5. Insert the folded Lung-Bag into the mouth of the manikin (Figure 5).
6. Push the Lung-Bag into the mouth until the bag exits through the neck of the manikin (Figure 6).
7. Grasp the Lung-Bag and finish pulling it through the neck of the manikin (Figure 7).
8. Attach the Face-Shield tabs to the 'buttons' located on the side of the manikin head (Figure 8).
9. Lightly pull down on the Lung-Bag to help open the mouth hole on the Face-Shield.
10. Then position the Lung-Bag on top of the compression plate (the hard plastic circular piece found on top of the torso foam core) with the circular imprint on the Lung-Bag aligning with the compression plate (Figure 9).
11. With proper positioning, the neck of the Face-Shield/Lung-Bag will have slack. Make sure it is not twisted.
12. Close the torso and secure the latch by pushing on the shoulder tabs until they snap into place.

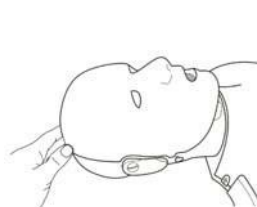


Figure 4

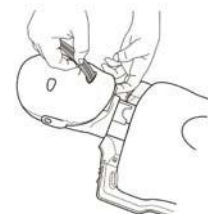


Figure 5

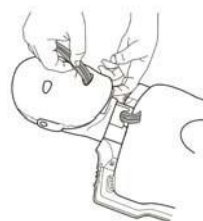


Figure 6

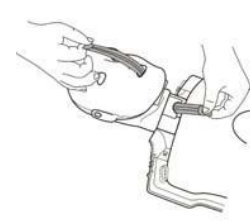


Figure 7



Figure 8

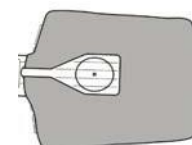


Figure 9

Prestan Professional ADULT MANIKINS



TROUBLESHOOTING GUIDE

SYMPTOM	POSSIBLE ISSUES OR COMMENTS	SOLUTIONS
CPR Rate Monitor		
Cannot find CPR Rate Monitor (Monitor) – no lights in the shoulder area	You have purchased a Prestan Professional Adult Manikin without the CPR Rate Monitor	If you were expecting a monitored manikin, check with the distributor where you purchased the manikin to insure you received what you paid for Note: CPR Rate Monitors may be purchased separately from your Prestan Authorized Distributor
CPR Rate Monitor (Monitor) lights are not working properly <i>To help understand how the Monitor works, please note that the Monitor lights are powered by a small current generated by the Clicker</i> <i>Please follow the troubleshooting tips at right in the order they are listed for best results</i>	- Chest compressions are not deep enough to signal Monitor - Batteries are not installed, need replacing or are not inserted correctly in the CPR Rate Monitor - The Clicker-Monitor connection is broken - Clicker may be clogged with plastic dust, interrupting the signal to the Monitor - Either the Clicker or the Monitor needs to be replaced	- Compress the chest until the click is audible from the Clicker at 2"-2.4" - Install or replace batteries into Monitor (back of manikin, left shoulder area) - requires 2-AA (1.5V) batteries - On the manikin's back, remove the long skinny Wire Cover between the Clicker and Monitor – squeeze connectors to disconnect the wires then reconnect to insure Clicker-Monitor communications - Remove Clicker and blow out dust – then replace Clicker and reconnect with Monitor - Exchange Clickers and Monitors with working manikins (if available) to determine whether the Clicker (RPP-ACLICK1) or the Monitor (RPP-AMON-1) needs to be replaced*
CPR Rate Monitor (Monitor) lights are all on	Batteries are running low	Replace batteries
CPR Rate Monitor (Monitor) red light is blinking	This is normal after replacing batteries or manikin is paused - may blink up to 16 times to indicate ready status	Commence or resume compressions - blinking will stop and feedback rate signals will begin

* Contact your Prestan Authorized Distributor to either purchase the replacement part or receive it at no charge if the manikin is within Prestan's 3-year warranty period and is a qualifying warranty repair.

Prestan Professional ADULT MANIKINS



TROUBLESHOOTING GUIDE (Continued)

SYMPTOM	POSSIBLE ISSUES OR COMMENTS	SOLUTIONS
Torso and Face-Shield/Lung-Bag		
Cannot hear a click when doing compressions	- Hand position and/or compression angle may be incorrect - Compressions may not be deep enough to activate the clicker - Clicker may be broken	- Confirm correct hand placement and compress straight down from directly above the manikin - Compress between 2 and 2.4 inches depth - Replace Clicker in back of manikin (RPP-ACLICK-1)*
Chest will not rise when breaths are given	- Head is not tilted back enough to open airway - Face Shield is not attached - Lung Bag may be damaged or torn (Note: there are holes in the Adult Lung Bag necks by design to allow deflation) - Lung Bag is not lying flat or properly positioned over the compression plate and/or neck of bag is twisted - Lung Bag has accumulated moisture in the neck, which may cause the sides to stick together and prevent air passage	- Push on forehead and/or lift at chin to tilt head back enough to open airway - Attach Face Shield at “buttons” on the side of the manikin’s head below ears - Blow into the Face-Shield/Lung-Bag while it is outside the manikin to make sure it inflates properly – if not, discard and use another Face-Shield/Lung-Bag - Revisit instructions for installing a Face-Shield/Lung-Bag earlier in this manual - A new Face-Shield/Lung-Bag should be installed after each class, or more often if needed due to excessive moisture buildup
Manikin torso will not stay closed Manikin torso opens when breaths are given	- Tabbed latches at shoulders may not be securely closed - Chest frame may be stretched from improper opening of torso - Tabbed latches at shoulders are broken from improper opening of torso	- Close tabs at shoulders – a snap will be heard when shoulder tabs are properly latched - Open torso and attempt to reverse stretched out frame by pushing it in the opposite direction - Torso Frame needs to be replaced (RPP-AFRAME-1)*

Prestan Professional ADULT MANIKINS



TROUBLESHOOTING GUIDE

(Continued)

SYMPTOM	POSSIBLE ISSUES OR COMMENTS	SOLUTIONS
Torso		
Rattling in chest	Compression Sleeve may be broken	Replace sleeve in chest with part RPP-ACOM-1 (compression assembly includes chest plate, spring, sleeve and clicker)*
Chest Plate is broken	Typically results from a combination of incorrect (not centered) hand placement, compression angle not straight down from directly above and/or overzealous compression force	Replace chest plate with RPP-ACOM-1 (compression assembly includes chest plate, spring, sleeve and clicker)*

* Contact your Prestan Authorized Distributor to either purchase the replacement part or receive it at no charge if the manikin is within Prestan's 3-year warranty period and is a qualifying warranty repair.

Prestan Professional ADULT MANIKINS



WARRANTY

Three Year Limited Warranty

1. **WARRANTY COVERAGE:** Prestan offers a LIMITED WARRANTY against any defects in workmanship or materials under normal use with respect to its PRESTAN TRAINING MANIKIN.
2. **WARRANTY PERIOD:** Prestan warrants the PRESTAN TRAINING MANIKIN for three (3) years from the date of purchase. Any component of the PRESTAN TRAINING MANIKIN which proves defective with thirty-six (36) months from the date of purchase will be repaired or replaced, at Prestan's option, free of charge to the owner. The owner is responsible for shipping and handling fees to the place which Prestan designates.
3. **WARRANTOR'S RESPONSIBILITIES:** Prestan will repair or replace (at Prestan's option) a defective warranted component subject to the provisions of Section 2 without charge except for transportation costs. OWNER'S SOLE AND EXCLUSIVE REMEDY AGAINST PRESTAN SHALL BE FOR THE REPAIR OR REPLACEMENT OF DEFECTIVE GOODS AS PROVIDED HEREIN. Repair or replacement of any component of any PRESTAN TRAINING MANIKIN does not extend this limited warranty or begin a new limited warranty period.
4. **CIRCUMSTANCE NEGATING WARRANTY AND OTHER LIMITATIONS:**
 - a. Prestan reserves the right to refuse warranty service or replacement when product failure is due to misuse, mishandling, misapplication, accident, neglect, improper installation, modification, use of unauthorized parts or attachments, adjustments or repair performed by anyone other than Prestan or its authorized agent, or any causes other than defective workmanship or material.
 - b. This limited warranty does not cover abrasion, abnormal use, damage caused by mishandling, neglect or unauthorized modification. This warranty is void if the PRESTAN TRAINING MANIKIN is structurally altered.
 - c. Owner must provide proof of the date of purchase, place of purchase and the purchase price in order to obtain service under this limited warranty.
 - d. Prestan makes no warranty other than as contained herein. No agent, representative, dealer or employee of Prestan has the authority to alter the obligations or limitations of this limited warranty.
5. **LIMITATIONS OF IMPLIED WARRANTIES AND DAMAGES:**
 - a. Any and all implied warranties on the PRESTAN TRAINING MANIKIN, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, SHALL NOT EXCEED IN DURATION THE RESPECTIVE TERMS OF THIS LIMITED WARRANTY WHICH BEGINS WITH THE DATE OF PURCHASE. Some states do not allow limitations on how long an implied warranty lasts, so the above limitation may not apply to you.
 - b. PRESTAN SHALL NOT BE LIABLE FOR INCIDENTAL OR CONSEQUENTIAL DAMAGES ARISING OUT OF THE USE OF THE PRESTAN TRAINING MANIKIN OR THE BREACH OF THIS OR ANY OTHER EXPRESSED OR IMPLIED WARRANTY. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation or exclusion may not apply to you. Any and all damages awarded for any breach of this limited warranty shall be limited to the amount of the purchase price.
6. **WARRANTY SERVICE.** To obtain warranty service contact your Prestan Authorized Distributor, or for information about Prestan, visit our website at www.prestanproducts.com.
7. **RIGHTS UNDER STATE LAW:** This warranty gives you specific legal rights and you may have other rights which vary from state to state.

Prestan[®]

AED UltraTrainer[™]



Prestan AED UltraTrainer™



Model PP-AEDUT-101 —> **Prestan AED UltraTrainer** (requires three (3) 'C' batteries; not included) preprogrammed with English/Spanish language, one set of adult/child training pads (pre-installed), one set of extra connector pins and a nylon carrying case.



Model PP-AEDUT-401 —> **Prestan AED UltraTrainer 4-Pack** (requires twelve (12) 'C' batteries; not included) preprogrammed with English/Spanish language, four sets of adult/child training pads, four sets of extra connector pins, 4 single nylon cases and a 4-pack nylon carrying case.



START HERE


- Use this instruction sheet as a reference for the Prestan AED UltraTrainer and Accessories.
- **Warning:** The Prestan AED UltraTrainer is intended for training purposes only. It is **NOT** an Automated External Defibrillator (AED) and will not deliver a shock. To avoid confusion, it should **NEVER** be kept in the same area with AEDs.
- Owner's information – For future reference, make note on your Quick Start Card of the Prestan Authorized Distributor where you purchased your AED Trainer and keep it with your unit. Contact your Prestan Authorized Distributor for product questions, reordering parts/accessories, or warranty issues.
- For information on Prestan's Authorized Distributors, help guides, resources and more, please visit the Prestan website at www.prestanproducts.com.

Features of the Prestan AED UltraTrainer


- High Performance Training Pads
 - ⇒ Clear Anatomical Placement Graphics
 - ⇒ Long-Life Adhesive
 - ⇒ Adheres well to all Manikins with no messy residue left behind
 - ⇒ Unique Pad Connector with locking Connector Pin
- Pad Sensing System
 - ⇒ Pads are pre-connected
 - ⇒ Pads on-manikin detection
 - ⇒ Remote Control NOT Required
 - ⇒ Plugging in connector NOT Required
 - ⇒ Works on any Manikin
- Protective & Reusable Pad Tray
- Adult / Child Selection Button
- Bright Illuminated Shock Button
 - ⇒ Flashes when shock is advised
- On/Off Button with Lighted Indicator and Program Pause Option
- Five Training Scenarios
- Clear & Calm Voice Prompts
- Chest Compression Rate Metronome
 - ⇒ Metronome can be turned on or off
 - ⇒ Metronome audible count reminders at 10, 20 and 30 compressions
- Automatic OR Semi-Automatic Shock-Delivery
- Programming Option to either include or not include time for breaths after compression cycle
- Long Battery Life
 - ⇒ Low Battery Voice Prompt Warning
 - ⇒ Automatic Shut Off Extends Battery Life
- Water & Splash Resistant
- Programming is compliant with current AHA Guidelines (2015)

ADULT/CHILD TRAINING PADS


- Pad cords are pre-connected to the AED Trainer
- Pad Sensing System automatically lets the unit know when the pads have been connected to a manikin
- Pads last for 100+ uses
- Unique pad replacement (see assembly section)
- Pads are stored on a convenient double-sided Pad Tray

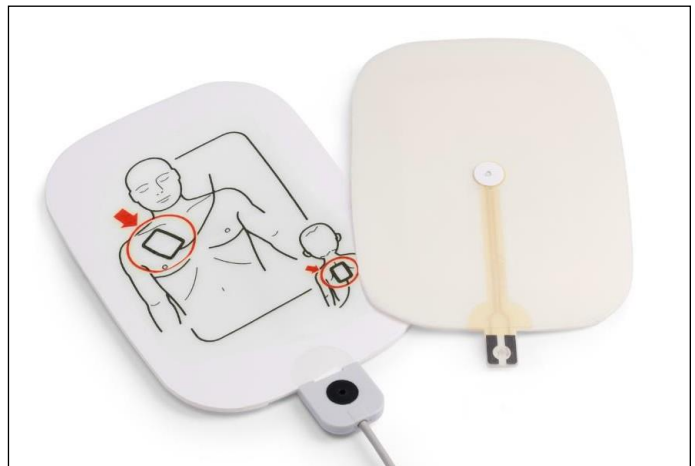

FEATURE GUIDE

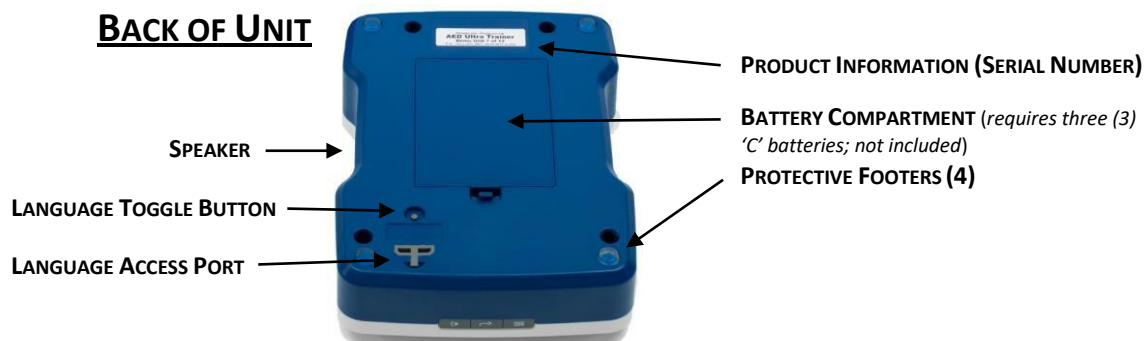

ACTION BUTTONS

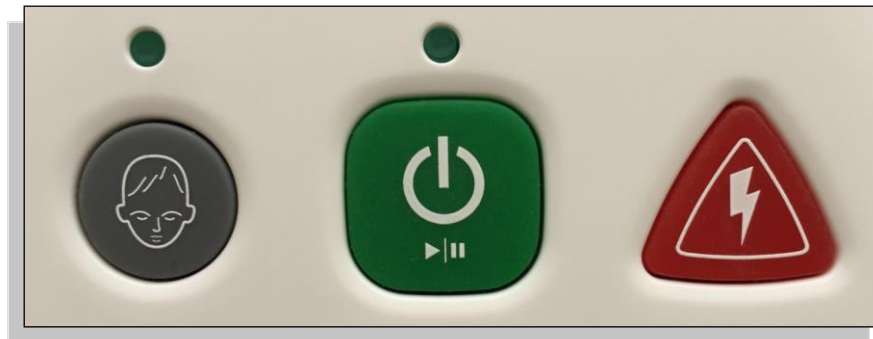

Child Button

On/Off/Pause Button

Shock Button

CHILD BUTTON

- After power is turned on, the programming offers an audible option to select the Child Button
- Indicator light above the Child Button will light up when the Child option has been selected and is active
- Indicator light turns off if the program is paused, but turns back on when the program is resumed

ON / OFF / PAUSE BUTTON

1. Push green button to start program / the indicator light above the button will turn on
2. Push green button again to pause program / indicator light will blink / push green button again to resume program
3. When unit is running, push and hold the button for 3 seconds to shut off unit

SHOCK BUTTON

- Red Shock Button will flash when program instructs to deliver shock
- In Semi-Automatic mode, pushing the Shock Button while it is flashing will stop the flashing and resume the program
- In Automatic mode, the Shock Button will flash approximately 6 times before auto-shocking and resuming the program.

OPTION BUTTONS



MENU BUTTON

The Menu Button allows customization of 4 options and can be used at any time. Push the Menu Button repeatedly to scroll through and hear each selection currently active for each of the options (see 1-4 below). Use the Selection Button to scroll through and change these options.

SELECTION BUTTON

Use with the Menu Button to change settings for each of the 4 Menu options. Repeatedly press Selection Button to hear current selection and scroll to the next. After you hear the option desired, simply stop pressing the button and that selection will remain active. A beep will be heard to confirm selection is set.

1. Available Training Scenarios:
 - a. Scenario 1 > Shock / Shock / Shock (default)
 - b. Scenario 2 > Shock / No Shock / No Shock
 - c. Scenario 3 > No Shock / No Shock / No Shock
 - d. Scenario 4 > Shock / No Shock / Shock
 - e. Scenario 5 > Shock / Shock / No Shock
2. Semi-Automatic / Automatic toggle
 - a. Semi-Automatic > Student must push flashing Shock Button to deliver shock (default)
 - b. Automatic > Unit will deliver shock when Shock Button is lit
3. Ventilations toggle
 - a. Ventilations On > Time is provided in scenarios to give breaths (default)
 - b. Ventilations Off > No time is provided in scenarios to give breaths
4. Metronome toggle
 - a. Metronome On > Audible ticks are heard for compression rhythm (default)
 - b. Metronome Off > No ticks are heard during compression section of program

See chart on Page 7 for Ventilations / Metronome interaction.

VOLUME BUTTON

- Toggle the Volume Button between 3 volume levels
- Works independently of Menu and Selections Buttons

FEATURE SETTINGS


	Ventilations On	Ventilations Off
CPR Metronome On	1. Program begins with analyzing patient and advising shock or no shock based on user scenario setting. 2. Audible prompt "BEGIN CPR" starts compression cycle. 3. Metronome beats and audible prompt notes "TEN", TWENTY, and THIRTY compression counts. 4. After 30 compression beats, the metronome pauses for an audible prompt to "GIVE 2 BREATHS". 5. Compressions beats and breaths continue for 5 cycles (2 minutes). 6. Unit analyzes patient. 7. Cycles repeat indefinitely following scenario shock pattern until unit is manually shut off.	1. Program begins with analyzing patient and advising shock or no shock based on user scenario setting. 2. Audible prompt "BEGIN CPR" starts compression cycle. 3. Compression beats with audible prompt to "CONTINUE CPR" every 24 seconds (x4). 4. Compression beats continue for 2 minutes. 5. Unit analyzes patient. 6. Cycles repeat indefinitely following scenario shock pattern until unit is manually shut off.
CPR Metronome Off	1. Program begins with analyzing patient and advising shock or no shock based on user scenario setting. 2. Audible prompt "BEGIN CPR" starts compression cycle. 3. Audible prompt notes "TEN", TWENTY, and THIRTY compression counts. 4. After time for 30 compressions, the metronome pauses for an audible prompt to "GIVE 2 BREATHS". 5. Compressions beats and breaths continue for 5 cycles (2 minutes). 6. Unit analyzes patient. 7. Cycles repeat indefinitely following scenario shock pattern until unit is manually shut off.	1. Program begins with analyzing patient and advising shock or no shock based on user scenario setting. 2. Audible prompt "BEGIN CPR" starts compression cycle. 3. Silence except for audible prompt to "CONTINUE CPR" every 24 seconds (x4). 4. Cycle stops at 2 minutes. 5. Unit analyzes patient. 6. Cycles repeat indefinitely following scenario shock pattern until unit is manually shut off.

CARE AND CLEANING INSTRUCTIONS




1. Store the Prestan AED UltraTrainer in its case and in a cool, dry area.
2. Return the training pads and cables to their respective storage pouches located on the interior lid of the unit's cover. If pads and cables are stacked on top of the unit, the cables may imprint grooves/indentations into the surface of the pads.
3. If the unit is being stored for long periods, remove batteries to protect against corrosion.
4. Never pick up or carry the unit by the cables.
5. The Prestan AED UltraTrainer should be kept clean. If it becomes soiled, simply wipe with a clean, damp, soft cloth. If the Trainer is extremely dirty, common cleaners can be used cautiously with a minimum amount of abrasion to ensure the life of the unit. Do not submerge unit in water.

PAD REPLACEMENT


1. Prestan Training Pads will generally last approximately 100 uses.
2. When pad replacement is needed, it is good practice to replace both pads at the same time.
3. Release the existing pads from the gray cable connectors by popping out the black pad connector pins with the tip of a ball point pen.
4. Slide the pads out from the cable connectors.
5. Insert new pads into the cable connectors.
6. Use the pad connector pins to lock the pads in place.

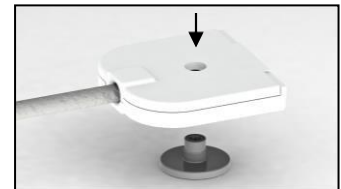

- Insert pad into cable connector (be sure the pad graphics are on same side as the large head of the connector pin and that the connector pin snugly fits into the indentation in the cable connector)

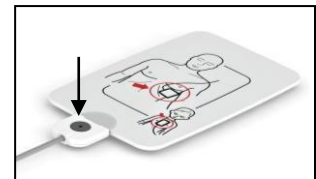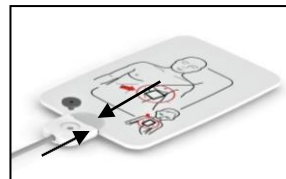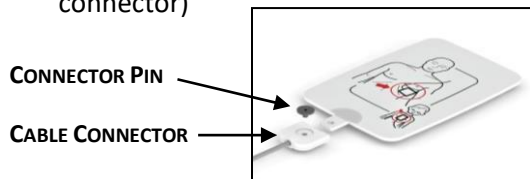

- Insure pad is not upside down
- If the pad has a protective backing over the adhesive (replacement pads), remove and dispose of this and adhere the pads to the double-sided pad case

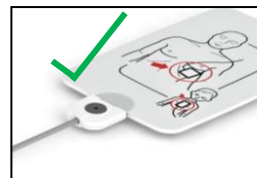

TROUBLESHOOTING GUIDE

SYMPTOM	POSSIBLE ISSUES OR COMMENTS	SOLUTIONS
Unit does not turn on	- Unit must be manually turned on - Batteries are not properly installed - Batteries are dead or are missing	- Push green power button - Re-orient batteries properly - Replace batteries
Both pads attached to manikin but device will not advance to analysis	- Operational solution - Pad is inserted upside down in the cable connector - Pad or Adhesive worn out	- Insure pads are pressed on firmly, specifically so that the sensor in the center of the adhesive is touching a surface - Push out black connector pin, flip pad around and re-secure with black connector pin - Replace pad (lasts 100+ uses)
Device proceeds to analyzing before both pads have been placed on the manikin	- Pads not properly positioned on tray - Look for pad riding up pad tray wall	Re-position pad on tray to lay completely flat so it is not touching side walls
Device delivers shock before shock button is pressed	Device is In Automatic mode	Change back to Semi-Automatic Mode
Device prompts "Replace Batteries"	Device will run 10 or more hours with low batteries	Replace batteries within 10 hours
Device suddenly shuts down or voice prompts are distorted	Batteries are very low, almost dead	Install new batteries
Pad pulls out from pad connector	Check for black connector pin missing	Install new black connector pin
Pad will not stick to the manikin	Adhesive is worn out	Replace pads (lasts 100+ uses)
Grooves/indentations are appearing on the surface of the pads	The cables are indenting into the surface of the pads from improper storage placement	Store cables and pads in their respective pouches in the interior lid of the unit's carry bag

WARRANTY

Three Year Limited Warranty

1. **WARRANTY COVERAGE.** Prestan (The Manufacturer) offers a LIMITED WARRANTY against any defects in workmanship or materials under normal use with respect to its PRESTAN AED ULTRATRAINER.
2. **WARRANTY PERIOD.** Prestan warrants the PRESTAN AED ULTRATRAINER for three (3) years from the date of purchase. Any component of the PRESTAN AED ULTRATRAINER which proves defective within thirty-six (36) months from the date of purchase will be repaired or replaced (at Prestan's option), free of charge to the owner. The owner is responsible for shipping and handling fees to the place which Prestan designates.
3. **WARRANTOR'S RESPONSIBILITIES.** Prestan will repair or replace (at Prestan's option) a defective warranted component subject to the provisions of Section 2 without charge except for transportation costs. OWNER'S SOLE AND EXCLUSIVE REMEDY AGAINST PRESTAN SHALL BE FOR THE REPAIR OR REPLACEMENT OF DEFECTIVE GOODS AS PROVIDED HEREIN. Repair or replacement of any component of any PRESTAN AED ULTRATRAINER does not extend this limited warranty or begin a new limited warranty period.
4. **CIRCUMSTANCE NEGATING WARRANTY AND OTHER LIMITATIONS.**
 - a. Prestan reserves the right to refuse warranty service or replacement when product failure is due to misuse, mishandling, misapplication, accident, neglect, improper installation, modification, use of unauthorized parts or attachments, adjustments or repair performed by anyone other than Prestan or its authorized agent, or any causes other than defective workmanship or material.
 - b. This limited warranty does not cover abrasion, abnormal use, damage caused by mishandling, neglect or unauthorized modification. This warranty is void if the PRESTAN AED ULTRATRAINER is structurally altered.
 - c. This warranty does not cover consumable items such as batteries that may be provided as a courtesy or training pads used beyond their indicated limit.
 - d. Owner must provide proof of the date of purchase, place of purchase and the purchase price in order to obtain service under this limited warranty.
 - e. Prestan makes no warranty other than as contained herein. No agent, representative, dealer or employee of Prestan has the authority to alter the obligations or limitations of this limited warranty.
5. **LIMITATIONS OF IMPLIED WARRANTIES AND DAMAGES**
 - a. Any and all implied warranties on the PRESTAN AED ULTRATRAINER, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE, SHALL NOT EXCEED IN DURATION THE RESPECTIVE TERMS OF THIS LIMITED WARRANTY WHICH BEGIN WITH THE DATE OF PURCHASE. Some states do not allow limitations on how long an implied warranty lasts, so the above limitation may not apply to you.
 - b. PRESTAN SHALL NOT BE LIABLE FOR INCIDENTAL OR CONSEQUENTIAL DAMAGES ARISING OUT OF THE USE OF THE PRESTAN AED ULTRATRAINER OR THE BREACH OF THIS OR ANY OTHER EXPRESSED OR IMPLIED WARRANTY. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation or exclusion may not apply to you. Any and all damages awarded for any breach of this limited warranty shall be limited to the amount of the purchase price.
6. **WARRANTY SERVICE.** To obtain warranty service contact your Prestan Authorized Distributor, or for information about Prestan, visit our website at www.prestanproducts.com.
7. **RIGHTS UNDER STATE LAW.** This warranty gives you specific legal rights and you may have other rights which vary from state to state.

PARTS & ACCESSORIES


Contact your Prestan Authorized Distributor for Pricing and Ordering




11678
Single Carry Bag for the AED UltraTrainer



11806
4-Pack Carry Bag for the AED UltraTrainer

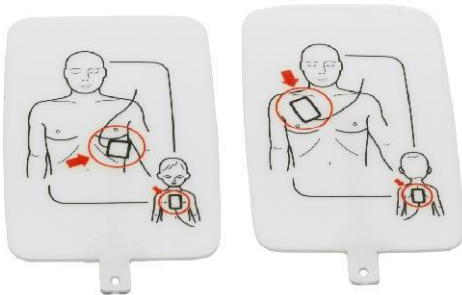

PP-UTPAD-1
Single set of Adult/Child Training Pads
for the AED UltraTrainer

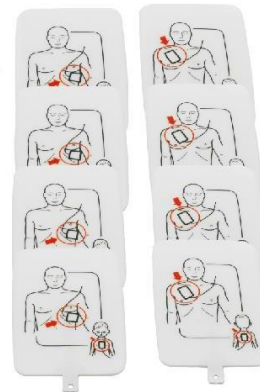

PP-UTPAD-4
Four (4) sets of Adult/Child Training Pads
for the AED UltraTrainer



RPP-AEDUT-CABLE
Replacement Cable Assembly for
the AED UltraTrainer



RPP-AEDUT-TRAY
Replacement Pad Storage Tray
for the AED UltraTrainer



PP-AEDUT-101-D English/Spanish
PP-AEDUT-102-D English/French
Language dongle for
the AED UltraTrainer

The Prestan AED UltraTrainer™

offers features not found on any other universal AED trainer.

This compact and lightweight unit gives students a realistic training experience with training pads that sense when they are placed on a manikin. In addition, the pads are preconnected and voice prompts are clear and calm.

Instructors can customize their training class with options to turn on or off the compression metronome and/or ventilation voice prompt.

With a 3-year warranty and durable long-lasting pads, the AED UltraTrainer is built to last.

Prestan is the First and Only Universal AED Trainer with:

- ✓ Dual Graphic (Adult/Child) Training Pads
- ✓ Pad Sensing System
- ✓ Preconnected Pads
- ✓ Child Training Button
- ✓ Child AED specific voice prompts
- ✓ Automatic or Semi-Automatic shock simulation
- ✓ CPR Coaching, "Give Breaths" voice prompt
- ✓ Turn on or off compression metronome
- ✓ Turn on or off "Give Breaths" prompts
- ✓ Low Battery warning
- ✓ 3 Year Warranty

Other features:

- ✓ Two languages (English, Spanish)
- ✓ 2015 Guidelines Compliant and Upgradable
- ✓ Play/Pause Button
- ✓ 5 Training Scenarios



Learn more about the Prestan AED UltraTrainer™
Scan the QR Code to watch the video

Prestan Products LLC © 2017



PRESTAN®

www.prestanproducts.com